

Thin Layer Chromatography (TLC)

Thin Layer Chromatography (TLC)


Introduction

- Thin Layer Chromatography (TLC) is a simple, rapid, and versatile chromatographic technique used for the separation and identification of components in a mixture.
- It is particularly important in pharmaceutical analysis due to its cost-effectiveness, speed, and reliability.
- TLC is a form of planar chromatography, where the stationary phase is supported on a flat surface (usually a glass or plastic plate).

Principle of TLC

- TLC works on the principle of adsorption and partition:
 - The stationary phase (e.g., silica gel or alumina) adsorbs the components of the sample.
 - The mobile phase (a suitable solvent or mixture) moves by capillary action across the plate.
 - The separation occurs due to differences in the affinity (adsorption and solubility) of each component toward the stationary and mobile phases.
- Components with higher affinity for the stationary phase move slowly, while those with greater solubility in the mobile phase move faster.

Technique of TLC

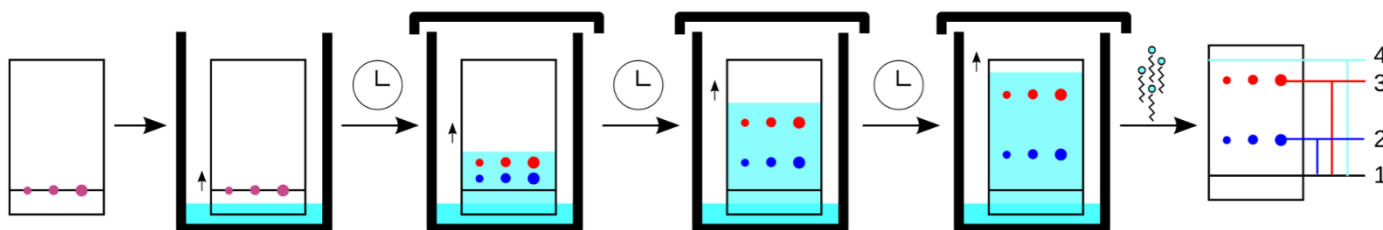

A. Preparation of the Plate

- A thin layer (0.1–0.25 mm) of adsorbent (silica gel, alumina, cellulose) is spread on a plate of glass, aluminum, or plastic.
- Plates can be commercially available pre-coated or manually prepared.

B. Activation

- The plates are heated in an oven at 100–120°C for about 30 minutes to remove moisture.

C. Sample Application

- Small spots (1–2 µL) of the sample solution are applied with a capillary tube or microsyringe near the bottom edge of the plate (~1 cm above the base).
- Multiple samples or standards can be spotted side by side.

D. Development of the Plate

- The plate is placed in a development chamber (e.g., a glass jar with a lid) containing the mobile phase solvent.
- The chamber is saturated with solvent vapor using a filter paper liner to ensure uniform movement.
- The solvent moves up the plate by capillary action, carrying the components with it.
- Development is stopped when the solvent front reaches near the top (~75% height of the plate).

E. Visualization

- After development, the solvent front is marked and the plate is dried.
- Spots are visualized using:
 - UV light (254 or 366 nm) if compounds are UV-active
 - Spray reagents (e.g., ninhydrin for amino acids, iodine vapors, sulfuric acid charring)
 - Fluorescence under UV if fluorescent indicator is present in the stationary phase

Retention Factor (Rf Value)

- The **Rf value** is a characteristic parameter used to identify compounds.

Formula:

$$R_f = \frac{\text{Distance traveled by the solute}}{\text{Distance traveled by the solvent front}}$$

Key Points:

- Rf is always between 0 and 1
- It is unitless
- It varies with solvent system, adsorbent type, thickness of layer, temperature, and humidity

Interpretation:

- Higher Rf: Indicates less polar compounds (move faster)
- Lower Rf: Indicates more polar compounds (stronger interaction with stationary phase)

Importance in Pharma:

- Comparison of Rf values of sample and reference under the same conditions helps identify unknowns.

Applications of TLC in Pharmaceutical Analysis

TLC is used both qualitatively and semi-quantitatively in pharma labs.

A. Identification of Drugs and Excipients

- Compare R_f of sample with that of reference standard.

B. Detection of Impurities and Degradation Products

- Useful in stability testing and impurity profiling.

C. Analysis of Herbal Drugs

- For fingerprint profiling of phytoconstituents.

D. Monitoring of Reactions

- Used in organic synthesis to track reaction progress.

E. Purity Testing

- Presence of additional spots indicates impurities.

F. Qualitative Estimation

- Estimation of content by spot intensity comparison with standard.

G. Bioanalysis

- Detection of metabolites in biological fluids (rare due to limited sensitivity).

H. Preparative TLC

- Scraping off and eluting separated components for further testing.

Advantages of TLC

- Simple and cost-effective
- Fast with no complex instrumentation
- Multiple samples can be analyzed simultaneously
- Useful for thermally unstable and non-volatile compounds

Limitations of TLC

- Low sensitivity compared to HPLC or GC
- R_f values may vary with experimental conditions
- Not suitable for quantitative analysis in routine pharma QC
- Manual method prone to operator error



FREE PREVIEW

[Click here to visit our website](#)

PHARMACY PREMIUM NOTES

Complete subject
and semesters are
available on
FirstHope app only.



[Scan or click here](#)

Click here to get the app



Available on Google Play

JOIN THE FIRSTHOPE COMMUNITY

Join our Telegram channel

Click here to join Telegram



[Scan or click here](#)

@FirstHopeNotes

Instagram

@firsthope.official

[Click here to follow >](#)

YouTube channel

FirstHope-Pharmacy / Live streams

[Click here to watch >](#)

WhatsApp channel

Follow FirstHope on WhatsApp

[Click here to join >](#)

FirstHope website

firsthopenotes.com

[Click here to visit >](#)

B.Pharm / D.Pharm / M.Pharm / Pharm.D

Preview notes. Complete subjects and semesters in the FirstHope app.