

Atomic Emission Spectroscopy

Atomic Emission Spectroscopy

- **Atomic Emission Spectroscopy (AES)** is an analytical technique used to determine the elemental composition of a sample by measuring the light (radiation) emitted by atoms when they return from an excited state to a lower energy state.

How it works:

1. The sample is atomized and excited using a high-energy source (like a flame, arc, spark, or plasma).
2. Excited atoms emit light at element-specific wavelengths as they relax back to their ground state.
3. The emitted light is dispersed and detected, and the intensity of the light is used to identify and quantify elements present.

Key Features:

- Used for qualitative and quantitative elemental analysis.
- Suitable for metal and non-metal detection in a variety of samples (solutions, solids, gases).
- Commonly applied in environmental analysis, metal alloys, biological samples, and industrial quality control.

Spectroscopic Sources in Atomic Emission Spectroscopy


Overview

- Spectroscopic sources are devices or methods used to excite atoms so they emit radiation (light).
- The emitted light is characteristic of each element, enabling both qualitative (what elements are present) and quantitative (how much of each element) analysis.

Purpose

- To provide sufficient energy so atoms in a sample transition from ground state to an excited state.
- When atoms return to ground state, they emit light at specific wavelengths, which are measured and analyzed.

Types of Spectroscopic Sources

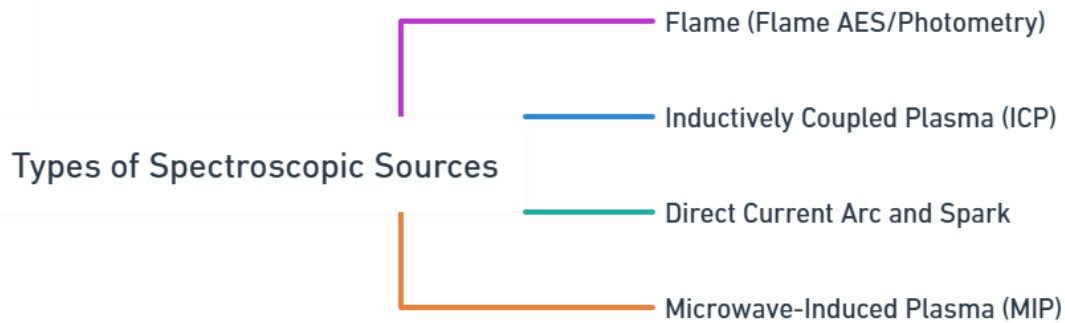

1. Flame (Flame AES / Flame Photometry)

- Uses a **fuel–oxidant flame** (e.g., air–acetylene, N₂O–acetylene) as a heat source.
- Best for **alkali and alkaline earth metals** (e.g., Na, K, Ca).
- Offers **limited excitation energy** → effective only for elements that are easily excited.

2. Inductively Coupled Plasma (ICP)

- Most widely used and powerful modern AES source.
- Employs an argon plasma generated by a radiofrequency (RF) induction coil (27–40 MHz).
- Plasma temperature: **~6000–10,000 K**.
- Excites a wide range of elements, including refractory and transition metals.
- **Advantages:**
 - High sensitivity
 - Multi-element detection
 - Minimal chemical interference

3. Direct Current (DC) Arc and Spark Sources

- Mainly for solid samples, especially metals.
- High-energy electrical discharge vaporizes and excites atoms directly.
- Less commonly used in pharmaceuticals due to lower precision and reproducibility compared to ICP.

4. Microwave-Induced Plasma (MIP)

- Plasma generated using microwave energy.
- Less common than ICP but useful in specialized applications, particularly when lower-cost plasma sources are needed.

Considerations


- Temperature stability of the source directly affects sensitivity and reproducibility.
- The source should minimize chemical interferences and ensure consistent excitation conditions.

Atomic Emission Spectrometers

Overview:

- Atomic Emission Spectrometers are instruments that detect and measure the intensity of light emitted by atoms excited by a spectroscopic source.

Main Components:

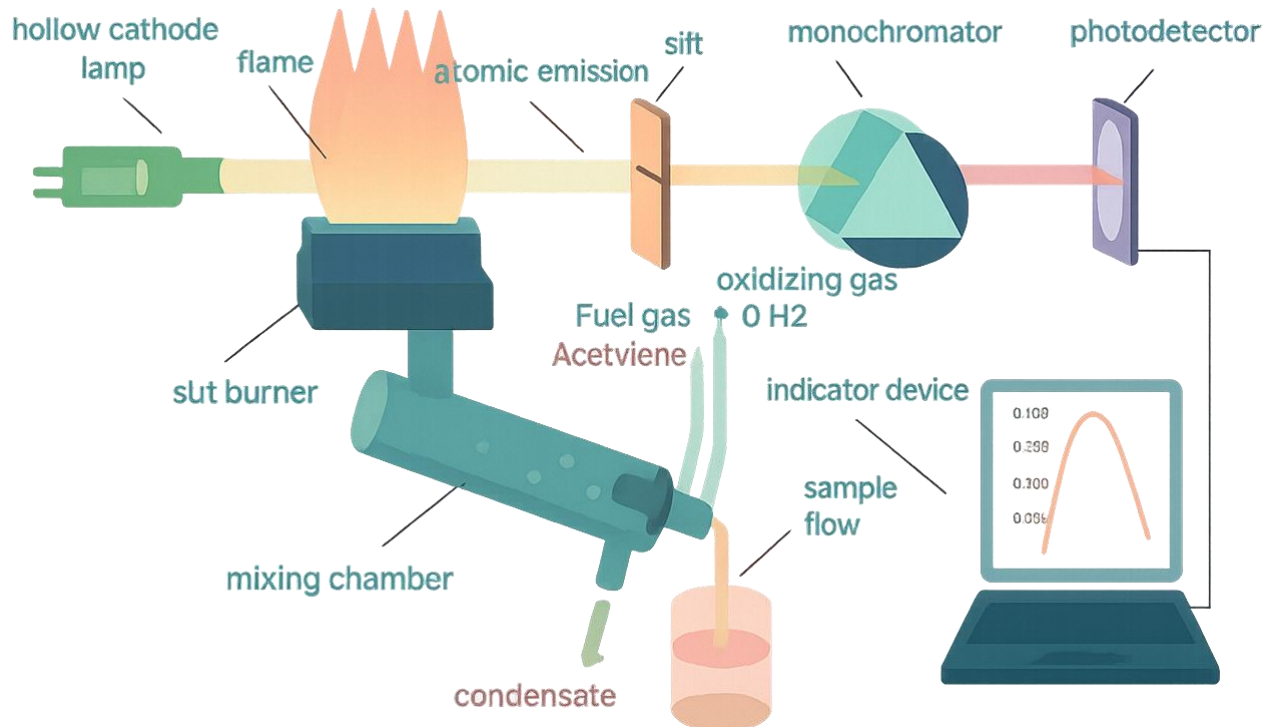

1. Sample Introduction System:

- Converts the sample (liquid or solid) into a form suitable for the excitation source.
- For example, a **nebulizer** in ICP converts liquid into aerosol.

2. Excitation Source:

- Flame, ICP, arc, or spark as discussed above.
- Excites atoms to emit characteristic radiation.

3. Wavelength Selector (Monochromator or Polychromator):

- **Monochromator:** Isolates a single emission line (wavelength) for detection.
- Uses prisms or diffraction gratings to separate light into its component wavelengths.
- Selects the most intense, interference-free line for each element.
- **Polychromator:** Allows simultaneous detection of multiple wavelengths for multi-element analysis.

4. Detector:

- Converts light into an electrical signal proportional to intensity.

- May be photographic (older methods) or photoelectric (modern methods).

5. Readout System:

- Displays or records the signal from the detector.
- Allows determination of **concentration** via calibration curves.

Performance Factors:

- **Resolution:** Ability to separate closely spaced emission lines.
- **Sensitivity:** Minimum detectable concentration.
- **Precision and Accuracy:** Reproducibility and closeness to true value.

Types of Spectrometers:



1. Sequential Spectrometers

- Measure one element at a time by scanning through different wavelengths.
- Good for limited number of analytes.

2. Simultaneous Spectrometers

- Measure multiple elements simultaneously.
- Use fixed-position detectors or photodiode arrays.
- Ideal for multi-element analysis in a single run (common in ICP-AES).

Advantages:

- High sensitivity and selectivity.
- Capable of multi-element analysis.
- Wide dynamic range.

Limitations:

- Some methods (e.g., flame AES) have limited sensitivity for elements that require higher excitation energies.

Photographic Detection in Atomic Emission Spectroscopy

Overview:

- One of the earliest detection techniques used in atomic emission spectroscopy.




FREE PREVIEW

[Click here to visit our website](#)

PHARMACY PREMIUM NOTES

Complete subject
and semesters are
available on
FirstHope app only.

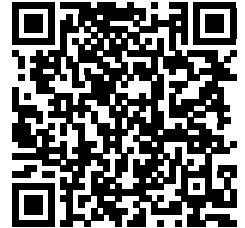

[Scan or click here](#)

Click here to get the app



Available on Google Play

JOIN THE FIRSTHOPE COMMUNITY

Join our Telegram channel

Click here to join Telegram



[Scan or click here](#)

@FirstHopeNotes

Instagram

@firsthope.official

[Click here to follow >](#)

YouTube channel

FirstHope-Pharmacy / Live streams

[Click here to watch >](#)

WhatsApp channel

Follow FirstHope on WhatsApp

[Click here to join >](#)

FirstHope website

firsthopenotes.com

[Click here to visit >](#)

B.Pharm / D.Pharm / M.Pharm / Pharm.D

Preview notes. Complete subjects and semesters in the FirstHope app.